
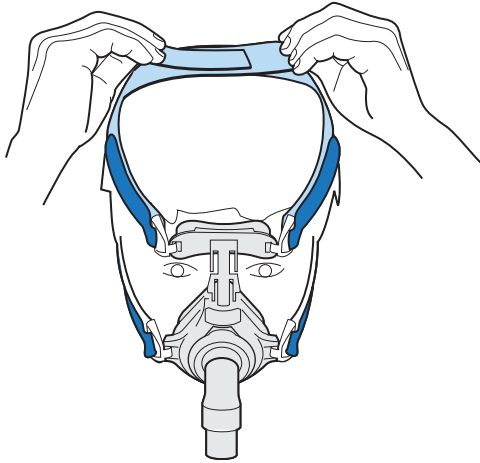
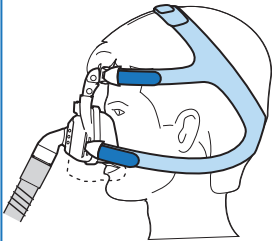
JOYCE SilkGel

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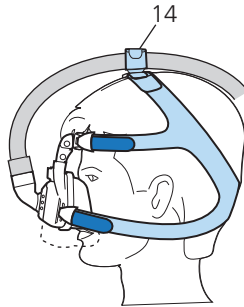Contents

| | | |
|-----------|----------------------------------|-----------|
| 1 | Introduction | 22 |
| 2 | Safety | 23 |
| 3 | Product description | 25 |
| 4 | Operation | 28 |
| 5 | Hygiene treatment | 30 |
| 6 | Disposal | 32 |
| 7 | Faults | 32 |
| 8 | Scope of supply | 33 |
| 9 | Technical data | 34 |
| 10 | Warranty | 37 |
| 11 | Declaration of conformity | 37 |

1 Introduction

1.1 Intended use

The JOYCE SilkGel nasal masks and full-face masks are used for treating sleep apnea and for the non-invasive ventilation of patients with respiratory insufficiency. They serve as a connecting element between the patient and the therapy device.

Use non-vented nasal masks only with an external exhalation valve.

Use non-vented full-face masks only in combination with therapy devices which have an active exhalation valve.

1.2 Contraindications

Erosions and ulcerations, skin allergies, skin rashes on the face, pressure points on the face, claustrophobia, anxiety, deformities of the face or nasopharynx, ingestion of drugs which may induce vomiting, the necessity for immediate intubation.

You should also observe the contraindications in the instructions for use for your therapy device.

1.3 Side effects

Blocked nose, dry nose, dry mouth in the morning, feeling of pressure in the sinuses, irritated mucous membrane, skin rashes, pressure points in the face, irritating noises when breathing.

2 Safety

Read these instructions for use through carefully. They are a constituent part of the mask and must be available at all times. In accordance with Directive 93/42/EEC, observe the points below.

2.1 Safety information

Risk of injury from the supply of oxygen!

Oxygen can become deposited in clothing, bed linen and hair. In conjunction with smoking, naked flame and electrical equipment, it can cause fires and explosions.

- ⇒ Do not smoke.
- ⇒ Avoid naked flame.
- ⇒ Do not use electrically-conductive tubes.

Risk of injury from oxygen therapy!

Oxygen therapies can lead to side effects if incorrectly used/ metered.

- ⇒ Only administer oxygen following prescription by a physician.
- ⇒ Follow the instructions for use of the oxygen system and the therapy device.

Risk of injury from re-inhaling CO₂!

If the mask is used incorrectly, CO₂ may be re-inhaled.

- ⇒ Do not seal the openings of the mask.
- ⇒ Only put mask on for a prolonged period if the therapy device is running.
- ⇒ Use non-vented nasal masks only with an external exhalation valve.

- ⇒ Use non-vented full-face masks only in combination with therapy devices which have an active exhalation valve.
- ⇒ Use the mask only within the quoted therapy pressure range.
- ⇒ Do not use mask on patients who cannot take off the mask independently or who tend to vomit.

Risk of injury from mask slipping!

If the mask slips or falls off, therapy is ineffective.

- ⇒ Monitor patients with restricted ability to breathe spontaneously.
- ⇒ Activate low pressure/leakage alarms on therapy device.
- ⇒ Subject patients with obstructive and restrictive lung diseases to individual monitoring.

2.2 General information

- If third-party items are used, functional failures and restricted fitness for use may result. Biocompatibility requirements may also not be met. Please note that in these cases, any claim under warranty and liability will be void if neither the accessories nor genuine replacement parts recommended in the instructions for use are used.
- The mask parts filled with gel are sensitive components. Handle them with care to prevent deformations and damage.
- To prevent an infection, bacterial contamination or functional impairments, follow the instructions in the section entitled "Hygiene treatment" (see 5, p. 30).

3 Product description

3.1 Overview

A diagram of the individual parts can be found on the enclosed folding sheet (**Figure A**).

Key

- 1 Headgear
- 2 Forehead cushion
- 3 Fine adjustment of the forehead support
- 4 Rough adjustment of the forehead support
- 5 Plug for pressure measurement connections
- 6 Pressure measurement connections
- 7 Retaining ring*
- 8 Elbow*
- 9 Rotating sleeve*
- 10 Emergency exhalation valve (on vented full-face masks only)
- 11 Mask body
- 12 Headgear clips
- 13 Mask cushion
- 14 Tube anchoring strap

* These parts are colored dark blue on non-vented variants (**Figure B**).

3.2 Compatible devices

You can use the mask with all therapy devices which are for the non-life-support treatment of patients with sleep apnea or respiratory insufficiency.

With many combinations of device, the actual pressure in the mask will not correspond to the prescribed therapy pressure, even if the therapy device is displaying the correct pressure. Have the device combination adjusted by a physician or specialist dealer so that the actual pressure in the mask corresponds to therapy pressure.

3.3 Exhalation system

Vented variants

The vented variants (which can be recognized by the transparent elbow and rotating sleeve) have an integrated exhalation system. The retaining ring and mask body are shaped so that a gap forms between these parts. The exhaled air can escape through this gap.

Non-vented variants

The non-vented variants (which can be recognized by the dark blue retaining ring, elbow and rotating sleeve, see **Figure B**) have no exhalation system.

Use non-vented nasal masks only with an external exhalation valve.

Use non-vented full-face masks only in combination with therapy devices which have an active exhalation valve.

3.4 Emergency exhalation valve

The vented full-face masks also have an emergency exhalation valve **10 (Figure C)**. The emergency exhalation valve has an opening to the atmosphere through which you can inhale ambient air if the therapy device fails. As long as the therapy device is in operation, the membrane is forced against the opening by the air flow and seals it. If the therapy device fails, the opening is free so that ambient air flows into the mask.

3.5 Pressure measurement connections

The mask has two connections **6** which you can use for measuring pressure or for supplying oxygen. If you are not using the connections, seal them using plug **5** in order to achieve therapy pressure (**Figure D**).

3.6 Quick-release feature (optional)

A quick-release feature is available as an option for full-face masks. The quick-release feature is for opening the headgear clips quickly in emergency situations. The quick-release feature consists of a special clip with a rip-cord which is fitted instead of one of the bottom headgear clips. When you pull on the rip-cord, the special clip opens immediately and you can remove the mask.

4 Operation

All the figures can be found on the enclosed folding sheet.

4.1 Put on/remove mask

1. Undo the bottom headgear clips on the mask.
2. Pull headgear **1** over your head and place the mask on the face.
3. Do up headgear clips **12**.
4. If the breathing tube is to be routed over the head: attach the tube to the headgear using tube anchoring strap **14** (**Figure I**).
5. Remove mask in reverse sequence.

4.2 Dismantle/assemble mask

1. Undo the headgear clips to remove headgear **1**.
2. If necessary: remove the plug for pressure measurement connections **5**.
3. Keep the bottom tabs of coarse adjustment **4** depressed and remove the forehead support (**Figure E**).
4. Keep the top tabs of coarse adjustment **4** depressed and remove fine adjustment **3** (**Figure E**).
5. Remove forehead cushion **2** from fine adjustment **3** (**Figure E**).
6. Take mask cushion **13** off mask body **11** (**Figure D**).
7. Turn retaining ring **7** counterclockwise and remove it (**Figure D**).
8. Take off elbow **8** through retaining ring **7**.
9. Take rotating sleeve **9** off elbow **8**.
10. If present: remove emergency exhalation valve **10** from the elbow (**Figure C**).

11. Assemble mask in reverse sequence.

You can tell whether the headgear is correctly attached by the hook-and-loop tabs. The hook-and-loop tabs on the forehead are light blue, those on the cheek are dark blue (**Figure F**).

4.3 Adjust mask (**Figure E**)

1. Keep the top tabs of coarse adjustment **4** depressed.
2. Push fine adjustment **3** into one of the three engaged positions.
3. If you need to adapt the height of forehead cushion **2**: rotate fine adjustment **3** through 180°.
4. Keep the bottom tabs of coarse adjustment **4** depressed.
5. Push coarse adjustment **4** into one of the three engaged positions (marked 1, 2 and 3) of the mask body.
6. Put on mask (see "4.1 Put on/remove mask" page28).
7. Adjust headgear **1** with the aid of the hook-and-loop tabs so that the mask seals on the face but is not too tight.

The hook-and-loop tab on top of the head is easiest to adjust by loosening the tab, gently pulling on the ends of the headgear with both hands and pushing the tab back down firmly (see **Figure G**).

5 Hygiene treatment

WARNING

Risk of injury if cleaning is inadequate!

Residues may block the mask, impair the integrated exhalation system and jeopardize the success of therapy.

- ⇒ For patients with a compromised immune system or particular background of illness, disinfect mask parts daily following consultation with the physician.
- ⇒ Thoroughly clean the contact surfaces between the elbow and the rotating sleeve.

5.1 Clean mask

1. Dismantle mask (see "4.2 Dismantle/assemble mask" page28).
2. Clean mask in accordance with the table below.

| Part | Daily | Weekly |
|--|----------|----------|
| Wash mask parts in warm water and a mild detergent. | X | |
| When washing mask parts (except for mask cushion and forehead cushion), clean thoroughly with a cloth or a soft brush. | | X |
| Wash headgear by hand. | | X |

3. Carefully rinse all parts with clean water.

NOTICE**Material damage from abrasion!**

Rubbing the mask cushion and forehead cushion can damage the coating.

⇒ Only clean mask cushion and forehead cushion when wet.

⇒ Do not rub mask cushion or forehead cushion when dry.

4. Allow all parts to air-dry.
5. Perform a visual inspection.
6. If necessary: replace damaged parts.

Discolored mask parts do not impair the function of the mask.

5.2 Tips

- You can wash all parts (except emergency exhalation valve and headgear) in a dishwasher at 65 °C.
- Do not dry parts in direct sunlight.
- Wash headgear before first use, as the color may run.
- Do not iron the headgear to keep the hook-and-loop tabs intact.
- Do not dry the headgear in a tumble-dryer.
- Handle the coating on the mask cushion and forehead support with care.

5.3 Change of patient

You must subject the mask to a hygiene treatment in the event of a change of patient. For information on a hygiene treatment in the event of a change of patient, see our brochure on the Internet (www.weinmann.de). If you do not have Internet access, contact Weinmann.

6 Disposal

You can dispose of all parts in domestic waste.

7 Faults

| Fault | Cause | Remedy |
|-------------------------------|--------------------|---|
| Pain in face due to pressure. | Mask is too tight. | Loosen headgear slightly. Modify position of forehead support. |
| Draught in the eyes. | Mask is too loose. | Tighten headgear slightly. Modify position of forehead support. |
| | Mask does not fit. | Use a different mask size. Contact your specialist dealer. |

| Fault | Cause | Remedy |
|----------------------------------|---|--|
| Therapy pressure is not reached. | Mask not adjusted correctly. | Readjust mask (see "4.3 Adjust mask (Figure E)" page29). |
| | Mask cushion damaged. | Replace mask cushion. |
| | Tube system leaking. | Check push-fit connector and fit of tube sleeves. |
| | Air escaping at pressure measurement connections. | Seal pressure measurement connections with plug (included in scope of supply). |

8 Scope of supply

| Size | Variant | Nasal mask | Full-face mask |
|-------------|-------------------|-------------------|-----------------------|
| S | vented | WM 26118 | WM 26418 |
| | non-vented | WM 26168 | WM 26468 |
| M | vented | WM 26128 | WM 26428 |
| | non-vented | WM 26178 | WM 26478 |
| L | vented | WM 26138 | WM 26438 |
| | non-vented | WM 26188 | WM 26488 |

9 Technical data

| | Nasal mask | Full-face mask |
|---|--|--|
| Device class as per Directive 93/42/EEC | IIa | |
| Dimensions (W x H x D) | approx. 80 mm x 130 mm x 90 mm ¹ | approx. 107 mm x 160 mm x 106 mm ¹ |
| Weight | approx. 85 g ¹ | approx. 102 g ¹ |
| Therapy pressure range | 4 hPa - 40 hPa | |
| Tube connection: tapered connection to EN ISO 5356-1 Vented variants: Non-vented variants: | Ø 22 mm (male) Ø 22 mm (female) | |
| Pressure measurement connections | Ø 4 mm | |
| Temperature range Operation Storage | +5 °C to +40 °C -20 °C to +70 °C | |
| Flow resistance Vented variants: at 50 l/min at 100 l/min Non-vented variants: at 50 l/min at 100 l/min Emergency exhalation valve: at 50 l/min at 100 l/min | 0.02 hPa 0.26 hPa 0.09 hPa 0.37 hPa -- -- | 0.02 hPa 0.11 hPa 0.09 hPa 0.37 hPa 0.5 hPa 2.0 hPa |

| | Nasal mask | Full-face mask |
|--|------------------------------|--------------------------|
| Quoted two-figure noise emission value to ISO 4871: - sound pressure level - sound power level (uncertainty factor: 3 dB(A)) | 15 dB(A) 23 dB(A) | 16.5 dB(A) 24.5 dB(A) |
| Switching pressure for emergency exhalation valve (vented variants) Open: Close: | -- -- | ≤ 1 hPa ≤ 2 hPa |
| Service life | up to 12 months ² | |

1 Depending on size and variant

2 The materials of the mask will age if they are exposed to aggressive detergents, for example. In individual cases, therefore, it may be necessary to replace mask parts sooner.

CE 0197

The right to make design modifications is reserved.


9.1 Materials

| Mask part | Material |
|--|---|
| Plug for pressure measurement connections | Silicone |
| Mask cushion, forehead cushion | Silicone and silicone gel, silicone coating |
| Mask body, elbow | Polycarbonate |
| Retaining ring, emergency exhalation valve | Polypropylene, TPE |

| Mask part | Material |
|-----------------------------------|---|
| Headgear clip, rip-cord clip | POM |
| Forehead support, rotating sleeve | Polyamide |
| Headgear | Lycra, polyester, polyurethane, UBL loop, cotton, polyamide |

9.2 Characteristic pressure/flow curve

The characteristic pressure/flow curve shows outlet flow rate as a function of therapy pressure.

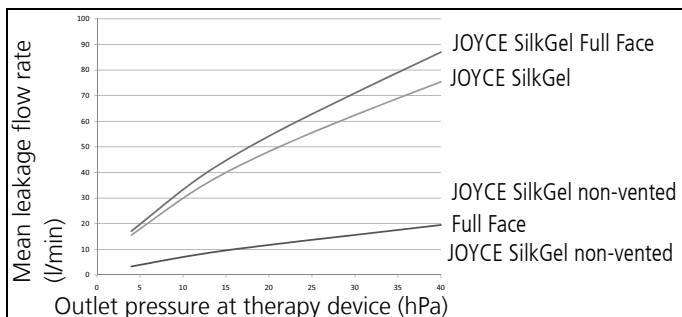10 Warranty

Weinmann gives the customer a limited manufacturer warranty on new genuine Weinmann products and any replacement part fitted by Weinmann in accordance with the warranty conditions applicable to the product in question and in accordance with the warranty periods from date of purchase as listed below. The warranty conditions can be downloaded from www.weinmann.de on the Internet. We can also send you the warranty conditions on request. In the event of a claim under warranty, contact your specialist dealer.

| Product | Warranty period |
|--|-----------------|
| Weinmann devices including accessories (except masks) for sleep diagnosis, sleep therapy, home ventilation, oxygen medicine and emergency medicine | 2 years |
| Masks including accessories, rechargeable battery, batteries*, sensors, tube systems | 6 months |
| Disposable products | None |

*Unless there is information to the contrary in the technical documentation.

11 Declaration of conformity

Weinmann Geräte für Medizin GmbH + Co. KG declares herewith that the product complies fully with the respective regulations of the Medical Device Directive 93/42/EEC. The unabridged text of the Declaration of Conformity can be found on our website at www.weinmann.de

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GmbH+Co. KG**